


Steerable Guiding Sheath

Guidee

Bi-Directional

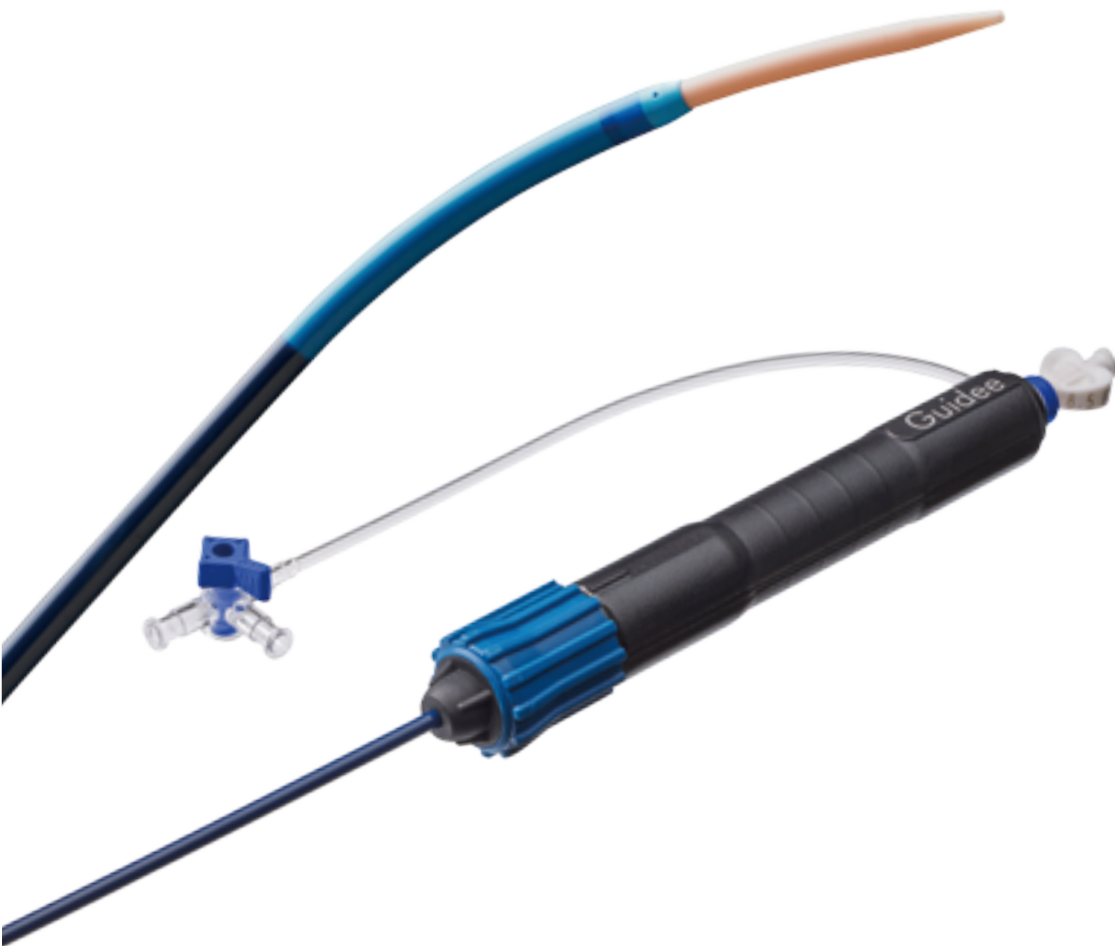

Guidee

Guide for Excellence

Precise Adjustment

The JLL designed dial structure provides fine maneuverability to reach the target site.

Neutral Position Marker

The neutral position marker helps physicians to find the neutral position visually and tactilely.

Smooth Transition

The tapered tip enables a smooth sheath and dilator transition, resulting in a gentle approach into the femoral vein and across the atrial septum with minimal resistance.

Soft Tip



The soft atraumatic tip reduces the risk of perforation.

Braided Shaft

The braided shaft enables a high flexibility and traceability.

Pre-Dilator

A 12F pre-dilation dilator comes with the sheath to provide an easier sheath insertion.

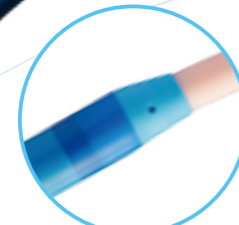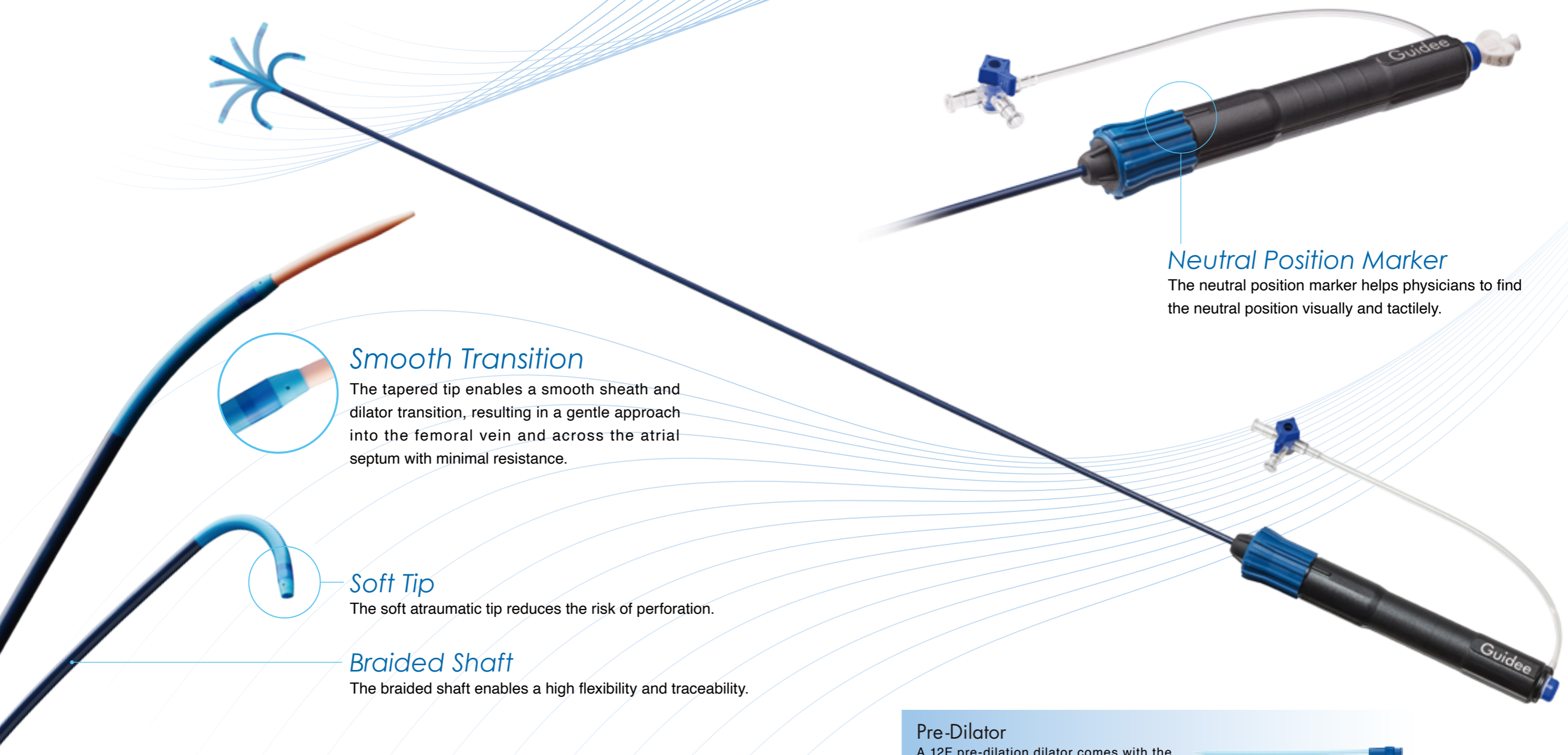

Specification

Steerable Guiding Sheath

Guidee

Model No.	Curve Diameter (mm)	ID (F)	OD (F)	Usable Length (cm)
ST8-64S	17	8.5	12	64

Accessories

50°Angled Dilator, Usable Length 85.3cm

0.032" J-Type Guidewire, Full Length 180 cm

12 Pre-Dilator, Usable Length 15.2cm

Some models may not be approved in your country. Please contact our local representatives for further information.

Japan Lifeline Co., Ltd.

2-2-20, Higashishinagawa, Shinagawa-ku,

Tokyo 140-0002 Japan

<https://www.japanlifeline.com>

2022-10-00-01-globalweb